

Can I output two different arbitrary waveforms with an SDG?


April 22, 2021


Yes.

Two-channel SDGs (SDG1000X, 2000X, and 6000X) can output two separate arbitrary waveforms.

1. Load the waveforms into the instrument
2. Select the channel you wish to load to and then press Save/Recall
3. Select the arbitrary waveform you wish to output for that channel
4. Repeat steps 2-3 for the other channel
5. Enable outputs

Here is an example of an SDG1032X with two separate arbs:

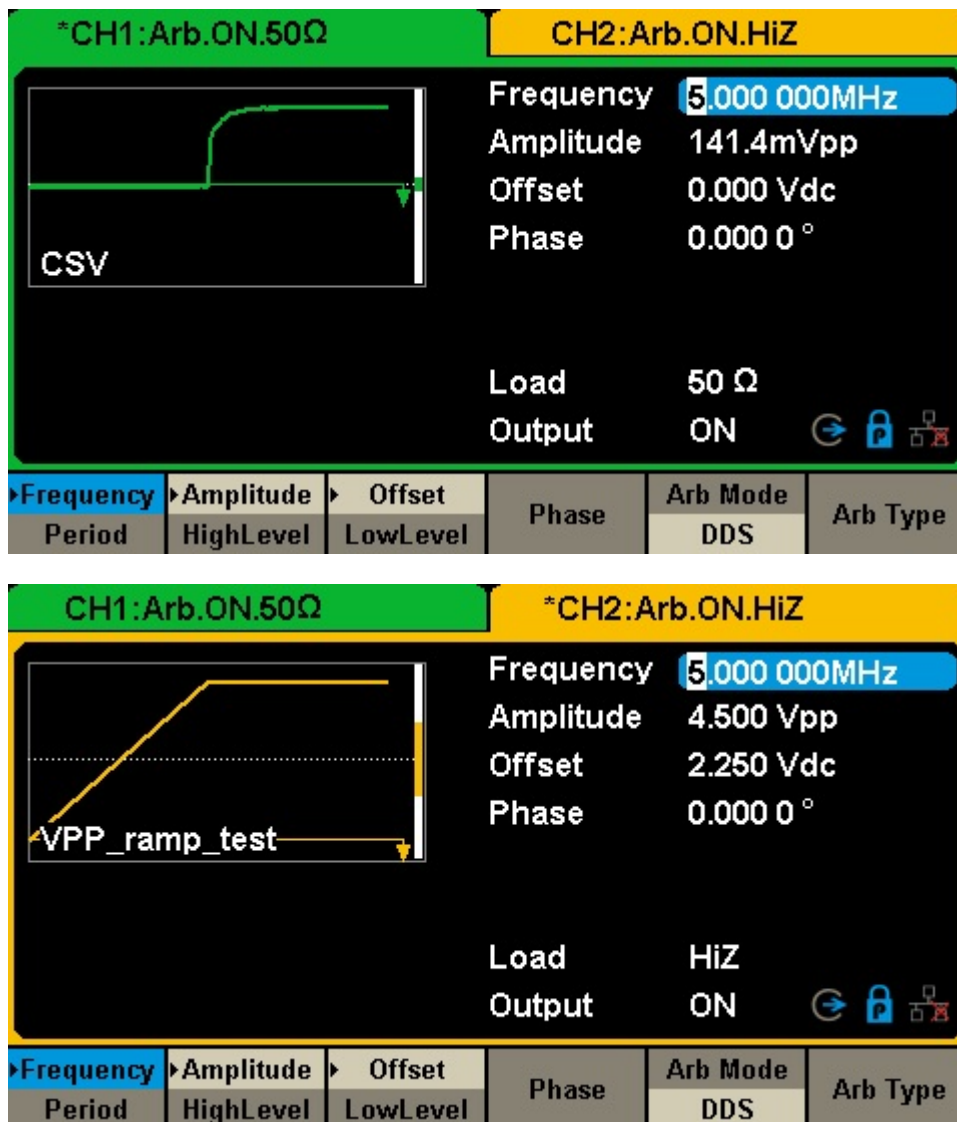

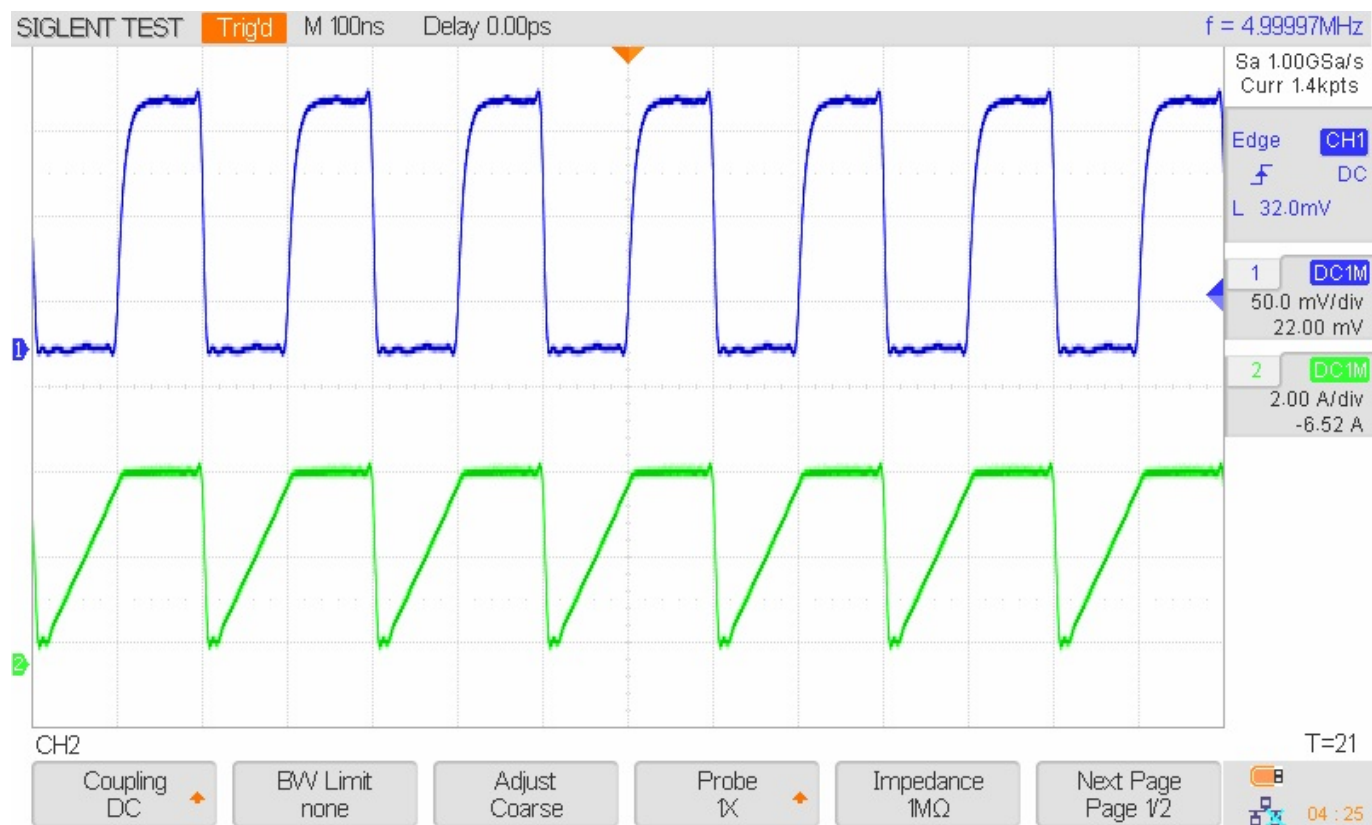



North American Headquarters

SIGLENT Technologies America, Inc
6557 Cochran Rd Solon, Ohio 44139
Tel: 440-398-5800
Toll Free: 877-515-5551
Fax: 440-399-1211
info@siglent.com
www.siglentamerica.com/

European Sales Offices

SIGLENT TECHNOLOGIES EUROPE GmbH
Staetzlinger Str. 70
86165 Augsburg, Germany
Tel: +49(0)-821-666 0 111 0
Fax: +49(0)-821-666 0 111 22
info-eu@siglent.com
www.siglenteu.com

Asian Headquarters

SIGLENT TECHNOLOGIES CO., LTD.
Blog No.4 & No.5, Antongda Industrial Zone,
3rd Liuxian Road, Bao'an District,
Shenzhen, 518101, China.
Tel: +86 755 3661 5186
Fax: +86 755 3359 1582
sales@siglent.com
www.siglent.com/ens