




How do I copy and move files from Local to a USB memory device with an analyzer?

September 10, 2020

1. Insert a FAT32 formatted USB stick into the front panel
2. Press FILE
3. Set Browser to DIR and use the knob to select Local. The local drive is the internal storage of the instrument. U-Disk0 is the USB device.

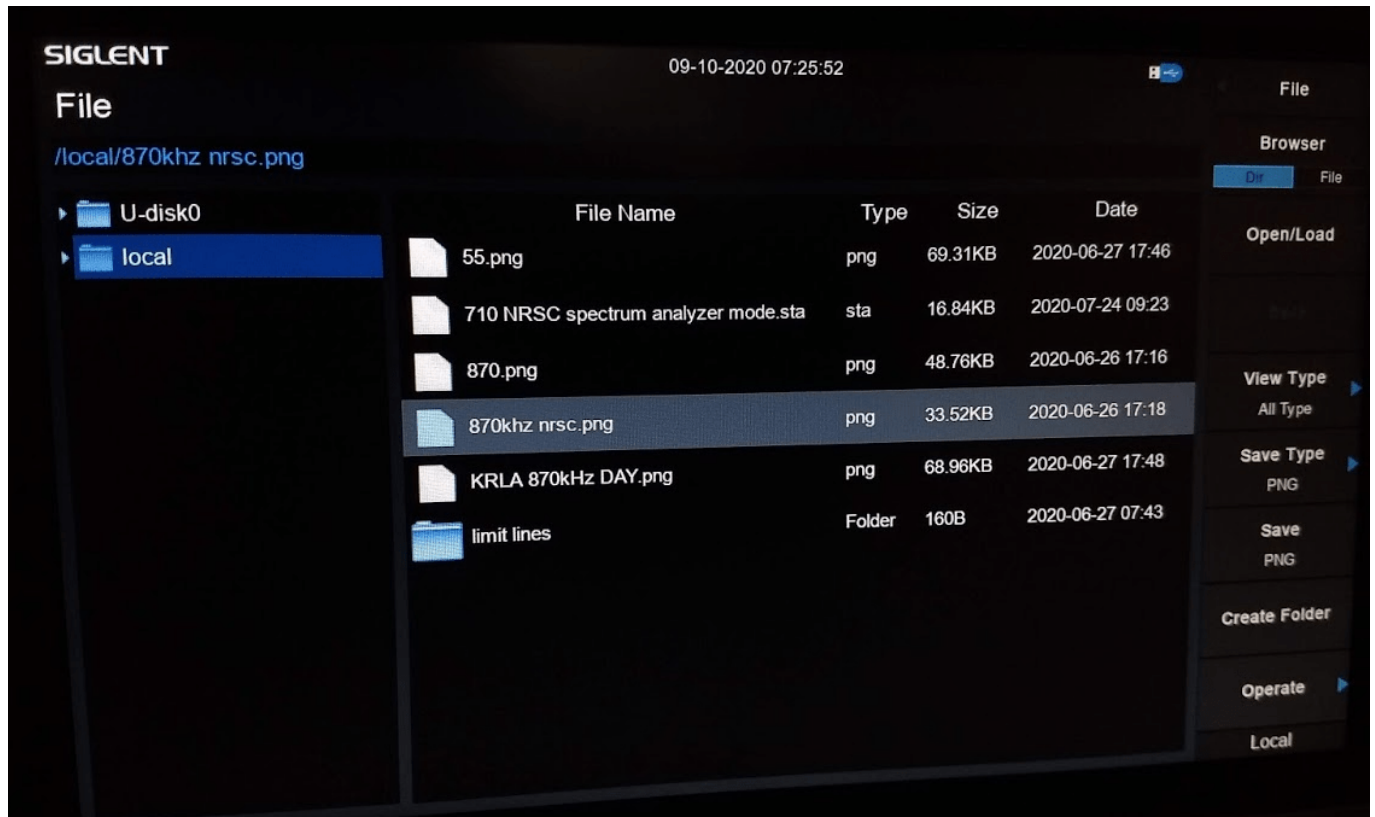


4. Use knob to highlight the file you wish to copy (or press Operate and "select all")

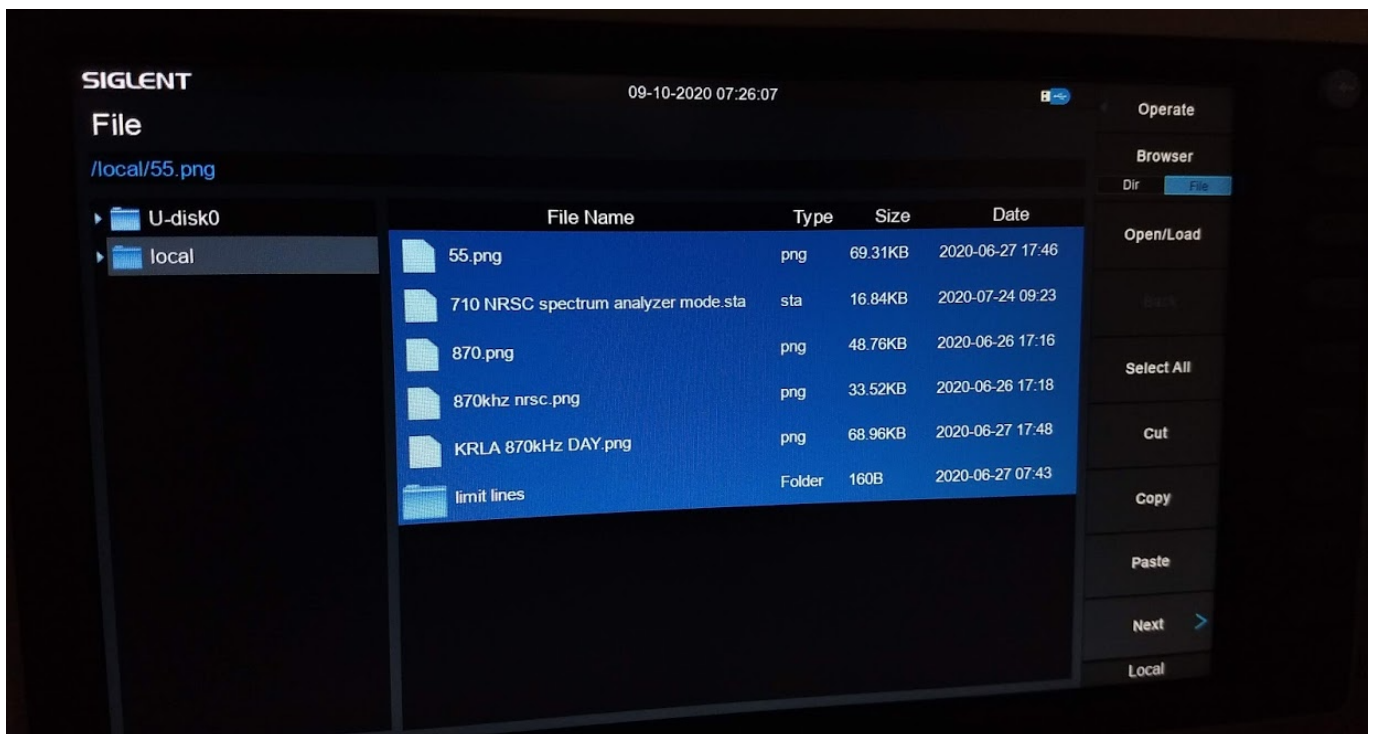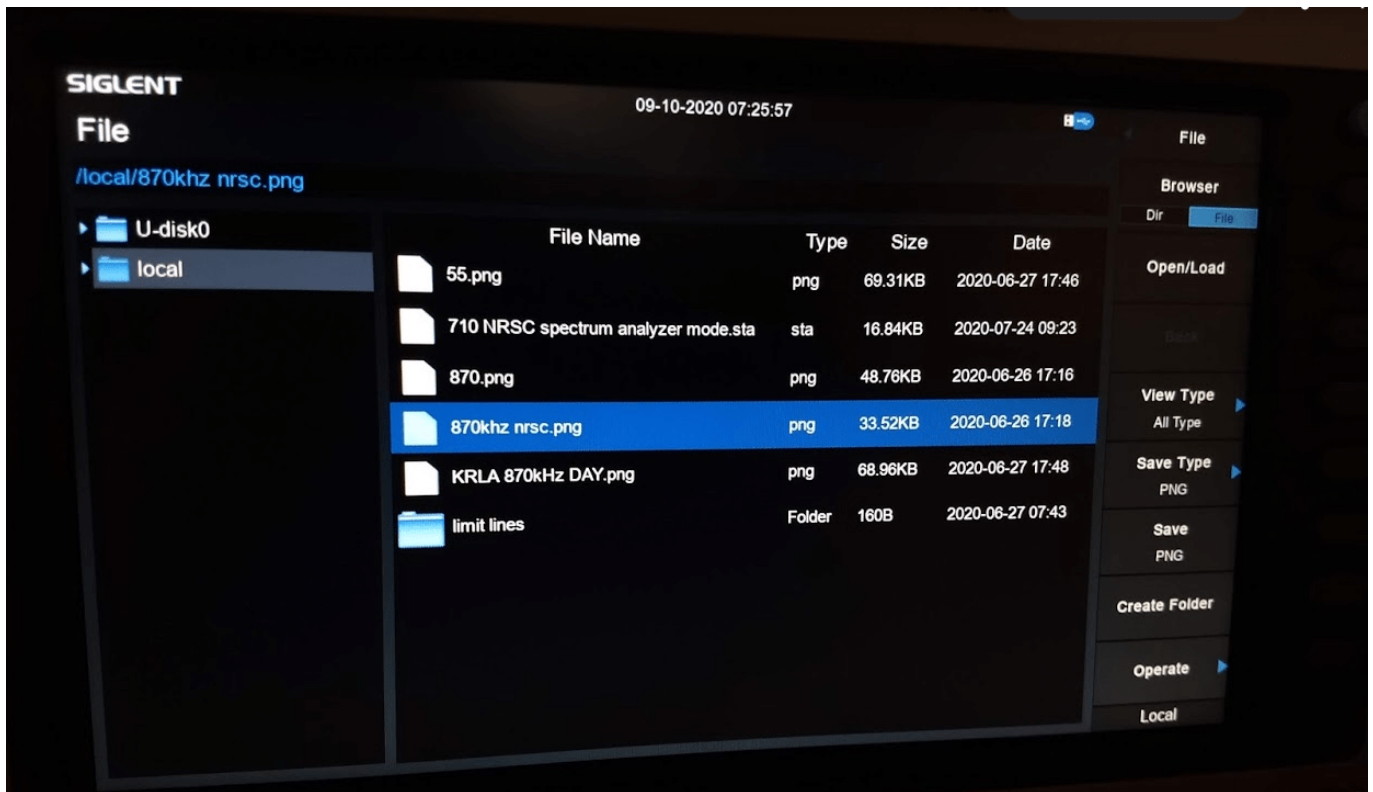


5. Press Operate > Copy

6. Press Browser > DIR and use the knob to select U-DISK

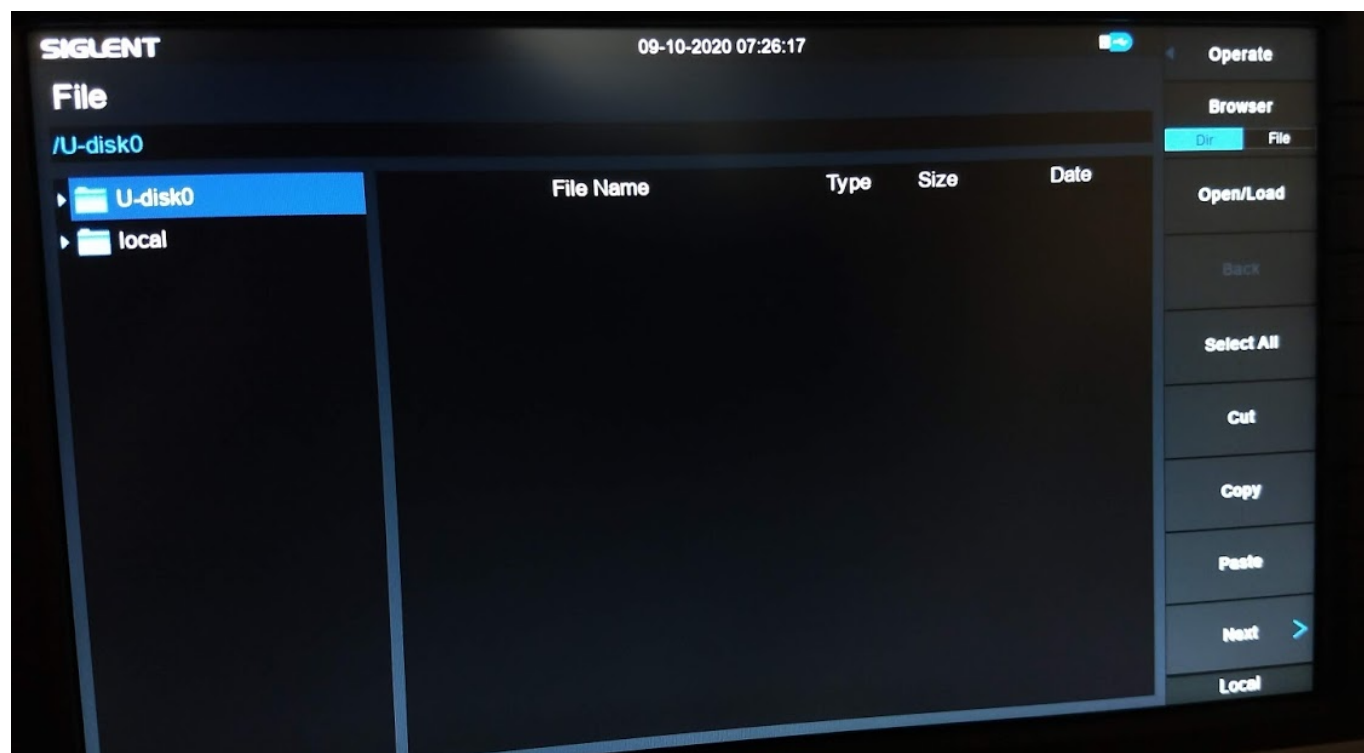

7. Press Paste

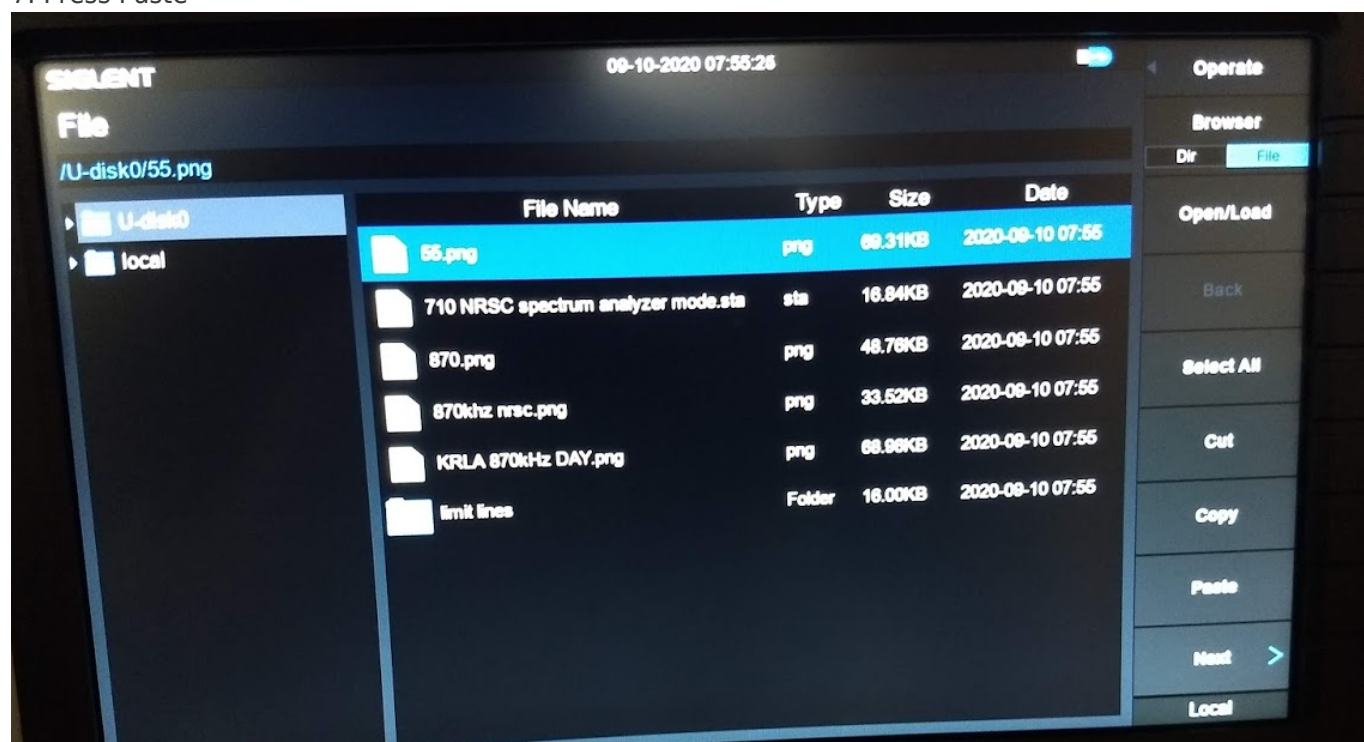



North American Headquarters

SIGLENT Technologies America, Inc
6557 Cochran Rd Solon, Ohio 44139
Tel: 440-398-5800
Toll Free: 877-515-5551
Fax: 440-399-1211
info@siglent.com
www.siglentamerica.com/

European Sales Offices

SIGLENT TECHNOLOGIES EUROPE GmbH
Staetzlinger Str. 70
86165 Augsburg, Germany
Tel: +49(0)-821-666 0 111 0
Fax: +49(0)-821-666 0 111 22
info-eu@siglent.com
www.siglenteu.com

Asian Headquarters

SIGLENT TECHNOLOGIES CO., LTD.
Blog No.4 & No.5, Antongda Industrial Zone,
3rd Liuxian Road, Bao'an District,
Shenzhen, 518101, China.
Tel: + 86 755 3661 5186
Fax: + 86 755 3359 1582
sales@siglent.com
www.siglent.com/ens